


BOOST

Technical sheet

English

Download
the Guide
at:

www.tartadesign.it

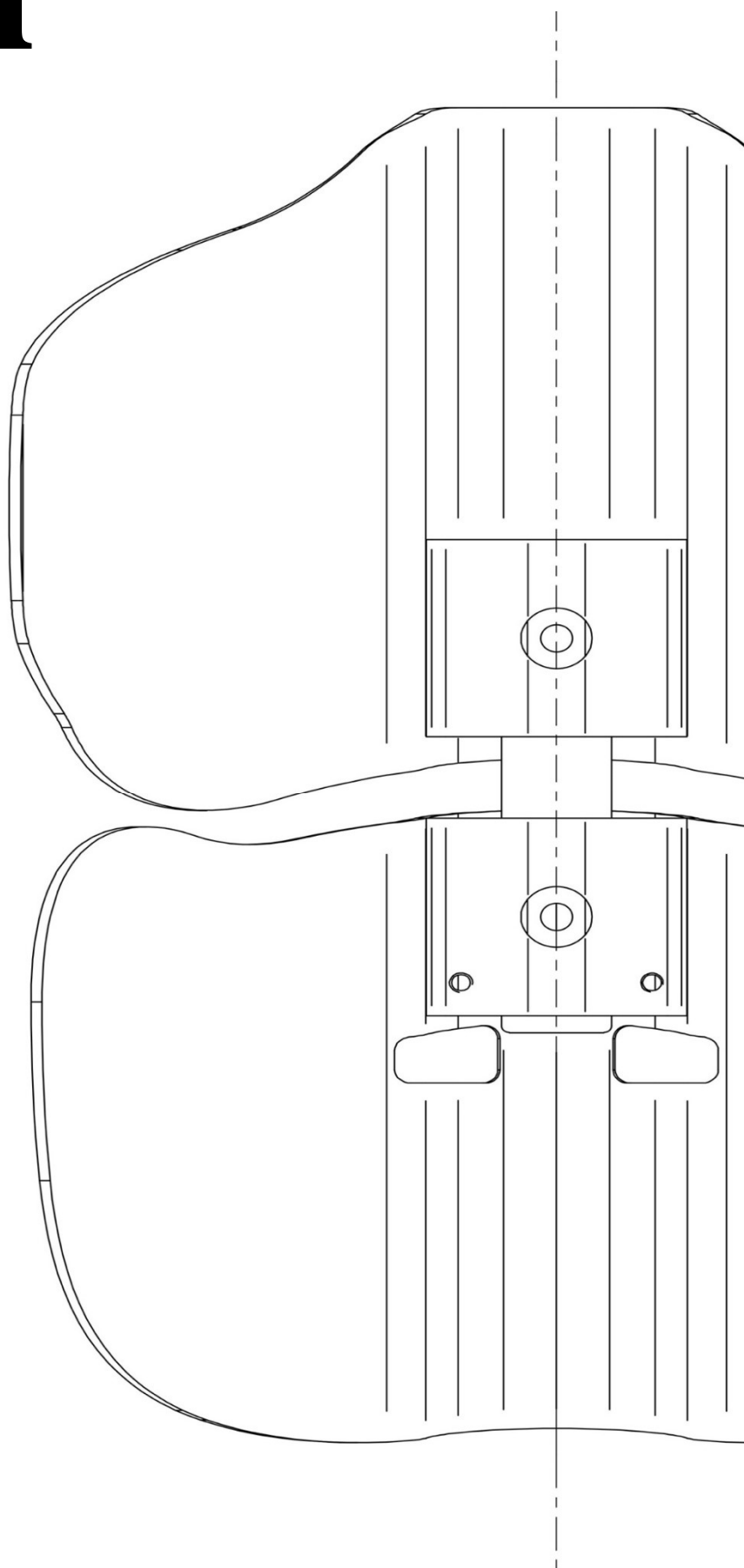

TARTA[®]
by Orthosanit Diffusion



Main features

Comfort

Each pad is made of foam, with a removable and washable covering produced with technical material for a maximum comfort.

Design

Minimal and soft lines highlight, an hi-tech design blending beauty and ergonomics.



Breathability

The space between the slats and the specific materials used for the pads, make the backrest highly breathable.

Flexibility

The central armonic spring steels and the modular system, allow the backrest to flex giving freedom of movement to the user.
New: with BOOST you can decide the level of flexibility of the backrest!

Height



▫ **LOW** ▫ **MID** ▫ **HIGH**

- **LOW** 27 cm
- **MID** 35 cm
- **HIGH** 43 cm

Size



▫ **M** ▫ **L** ▫ **XL** ▫

- **M** for wheelchairs from 35 to 40 cm
- **L** for wheelchairs from 38 to 43 cm
- **XL** for wheelchairs from 40 to 45 cm

Depth



STANDARD

5,5 cm



CONTOUR







9 cm

Aluminum is the most adaptable material to make each slide of the backrest suit the user's body!





Models

	STANDARD	CONTOUR
LOW M-L-XL		
MID M-L-XL		
HIGH M-L-XL		

How to bend the staves?

1 Press



2 Press to bend the stave. Maximum excursion 3 cm.


3 heights
3 sizes
2 depths

18 models










Adjustable clamp for rigid frame wheelchairs

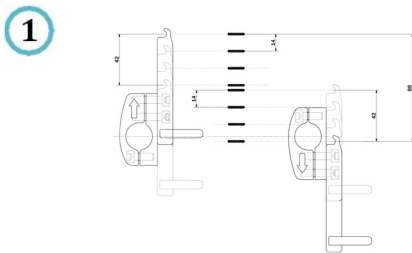





Features:

- 1 height adjustment
- 2 rotation adjustment
- 3 depth adjustment
- 4 tube diameter: 19-22-25mm
- 5 easy unhook
- 6 safety rubber band

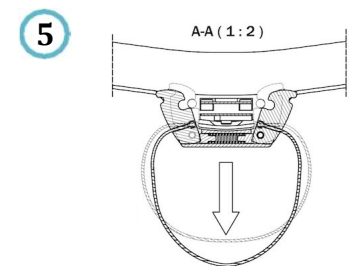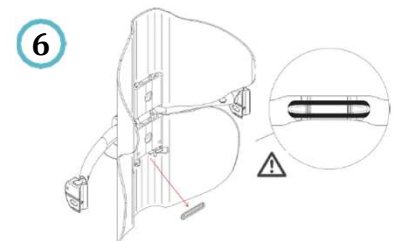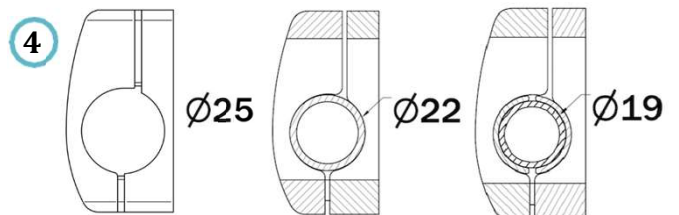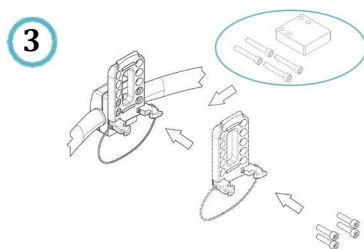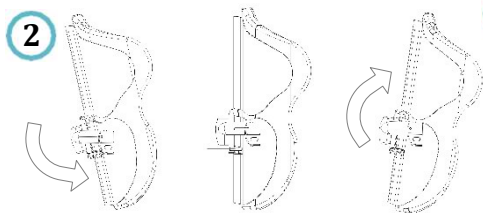








Adjustable clamp for folding wheelchairs

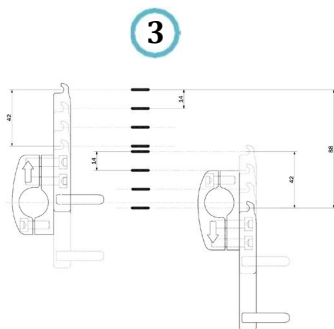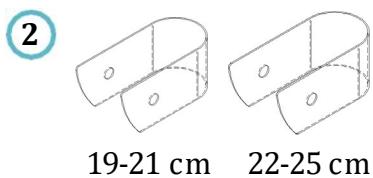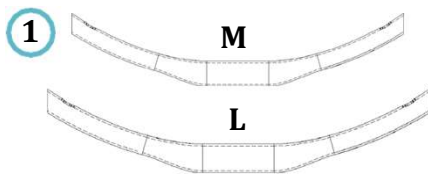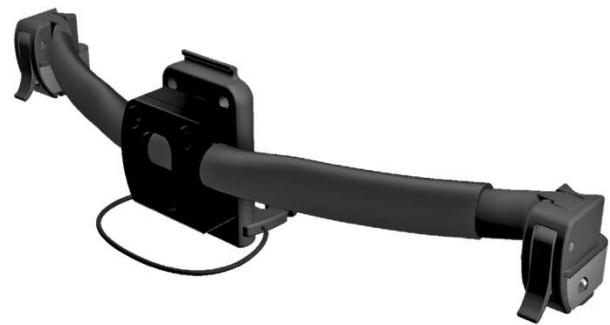

Features:

- ① rigidizer bar to specify:
M: 36 - 41 cm L: 41 - 46 cm
- ② hose clamp: 19-21 / 22-25 mm
- ③ height adjustment
- ④ depth adjustment
- ⑤ rotation adjustment
- ⑥ easy unhook
- ⑦ safety rubber band

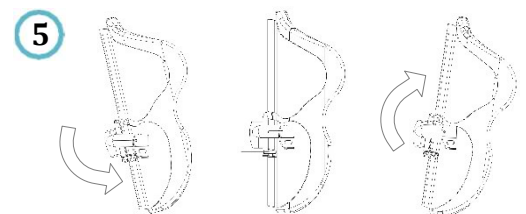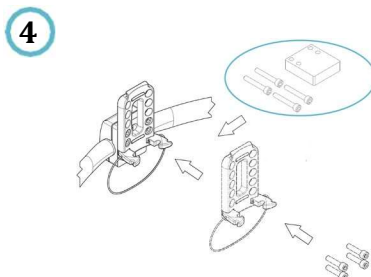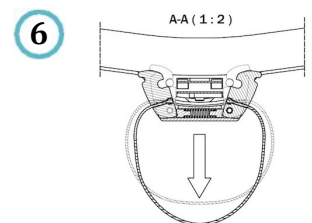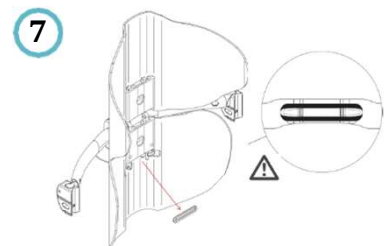

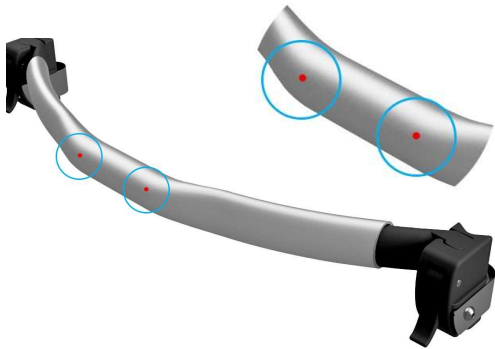

Balance points

on the rigidizer bar to simplify the balance of the clamp for folding wheelchairs.

2 rubber tubes (internal diameter 22mm)
easy adaptable to every wheelchair **to protect the backrest.**




Highly breathable pads

removable and **washable** at 30°.



Choose your **style.**



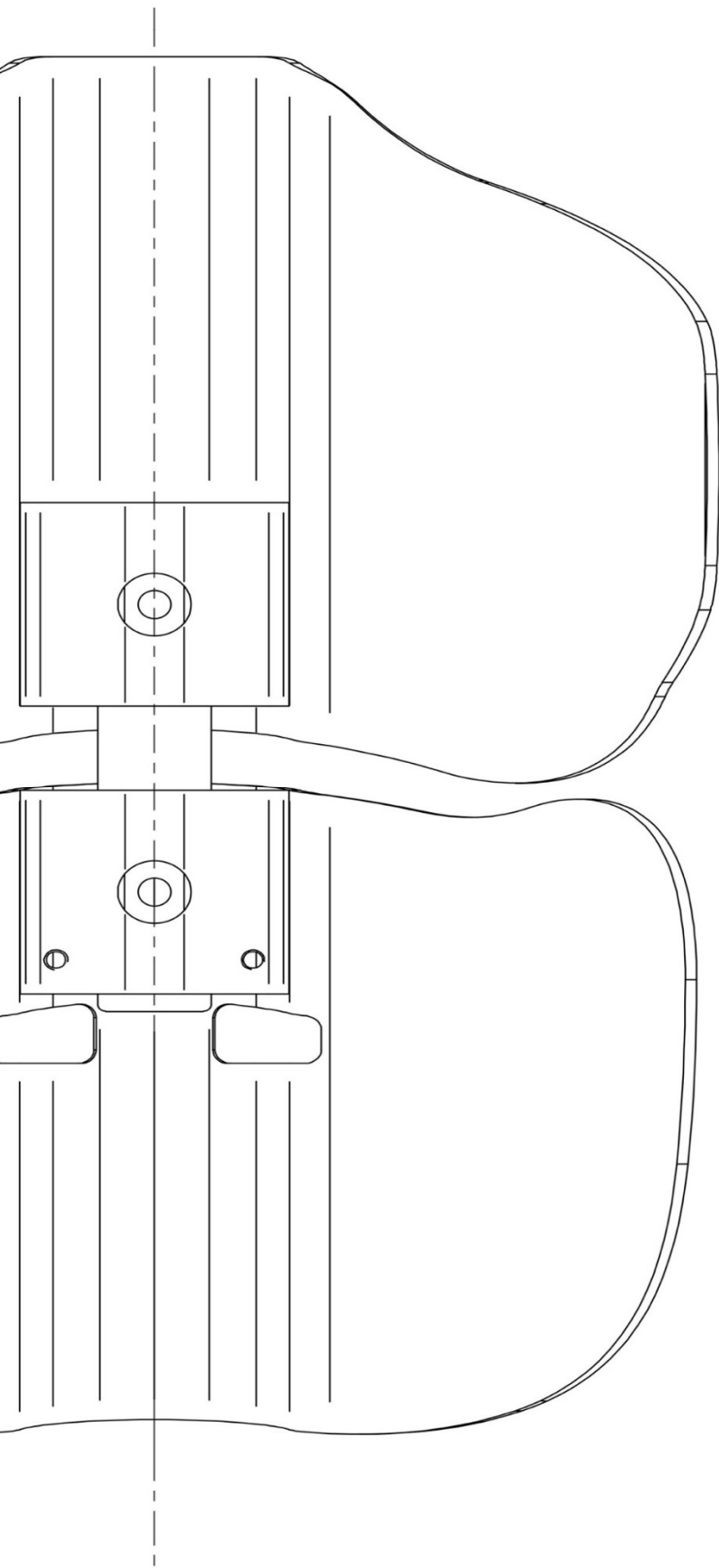

Questions?
Info?
Help?

Visit our website:
www.tartadesign.it

or contact us at:
info@tartadesign.it

Tarta by Orthosanit Diffusion

Registered Office:

Galleria Goito, 18, 19121

La Spezia, Italy

Headquarter:

Via Cotonificio, 56, Pesian di

Prato, 33037 (UD) Italy

www.tartadesign.it

info@tartadesign.it